



Superior Fluid Lock-In with Gelling Fibre Technology


Improved Wound Care Management



Kotszr™

Simplifying Wound Care

Kotszr™ Carboxymethyl Cellulose Dressings

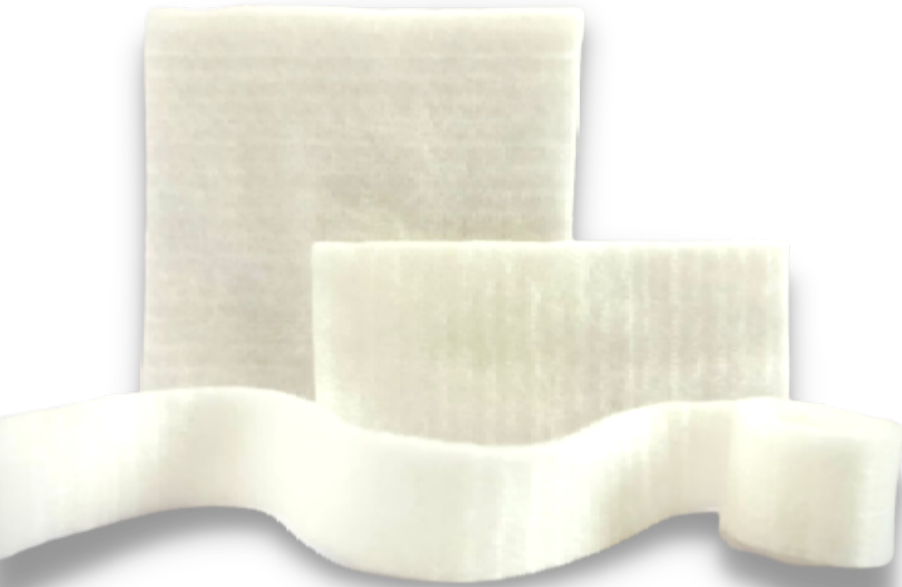


- *Kotszr™ CMC*, designed to provide **High absorbency** and **Tensile Strength** for **enhanced** patient comfort.
- Technologically advanced blend of **Carboxymethyl cellulose** and **Non-modified cellulose fibres**, to ensure accelerated wound healing.
- *Kotszr™ CMC* the dressing has been designed with a **unique construction** that allows for **intact removal** and **easier** dressing changes.


Key Features




Fluid Lock-In

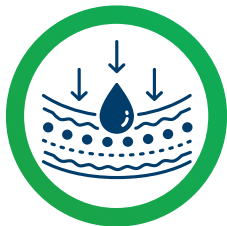

**High Absorbency with
Excellent Tensile Strength**

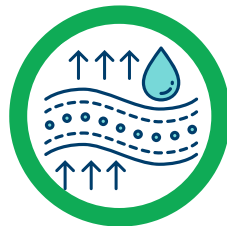

**High Vertical Absorption
Limits the Spread of Fluid
to Wound Edges**

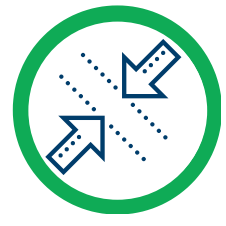

**Can be used
Under Compression**




**Protects periwound area
from Maceration**




Soft and Conformable

Mechanism of Action

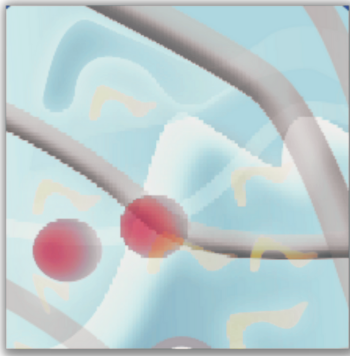


Fluid Lock-In

- **Traps** wound exudate, microorganisms & **reduces lateral spread** of fluid.
- High exudate absorption, quick **gel formation** by **exchange of ions**, helps to **protect wound** from maceration.

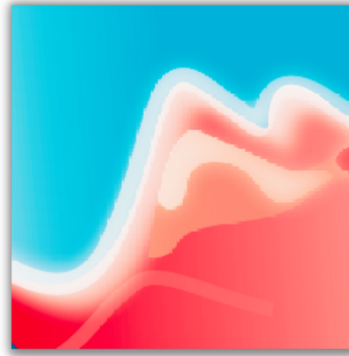

Micro Contours


- **Moulds** to the wound bed leaving **no dead space** for microorganisms to grow.
- **Holds debris** and **microbes** contained in exudate, ensuring a **clean wound**.

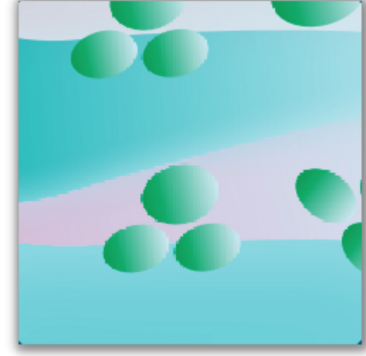

Responds

- Alters as per wound Fluid/Exudate levels by forming a **cohesive gel**.
- **Aids in wound healing** as microbes remain in the dressing **preventing colonization**.

Indications

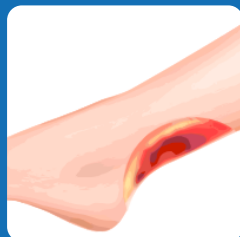


Venous Leg Ulcer



Burn



Pressure Ulcer

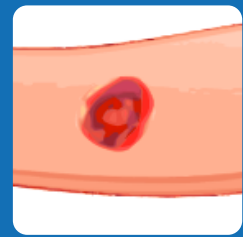


Cavity Wound



Donor Site



Diabetic Foot Ulcer

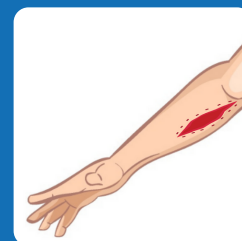

Abrasion




Laceration

- Superficial or Deep Wounds
- Moderately to Heavily Exuding Wounds
- Traumatic Wounds (e.g. Abrasions and Lacerations)

In Vitro Data Analysis


ABSORBENCY

- *Kotszr™ CMC* dressings have demonstrated **greater absorbency** compared to other commercially available brands.
- **Good Exudate Management** promotes positive healing progression of the wound.

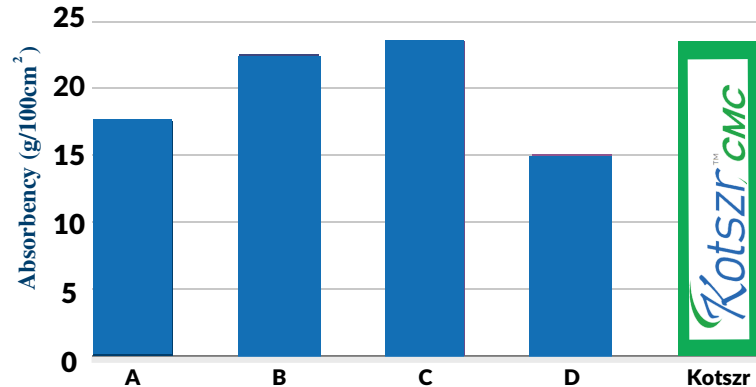


WET TENSILE STRENGTH

- *Kotszr™ CMC* dressings have significantly **higher wet tensile strength** when compared to other commercially available brands.
- It aids in **easy, trauma-free removal** of dressing from the wound site.

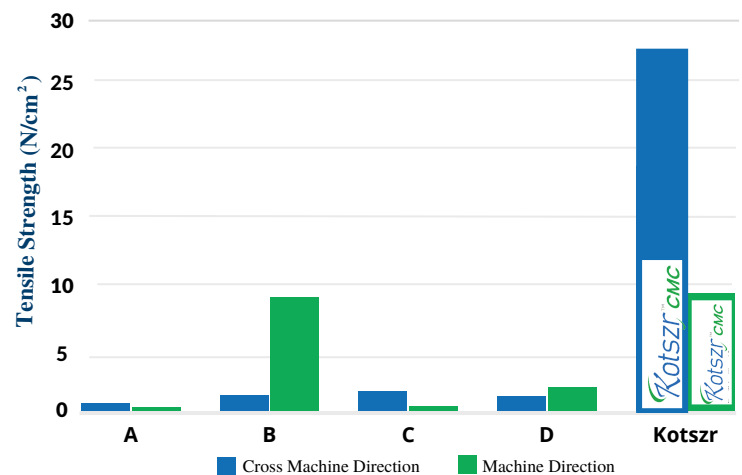

LATERAL WICKING (MM)

- *Kotszr™ CMC* dressings have the ability to **immobilize** absorbed fluids and prevent the **lateral spread** of exudate.
- *Kotszr™ CMC* dressing **controls fluid circulation** in the wound, **avoiding peri-wound maceration**.

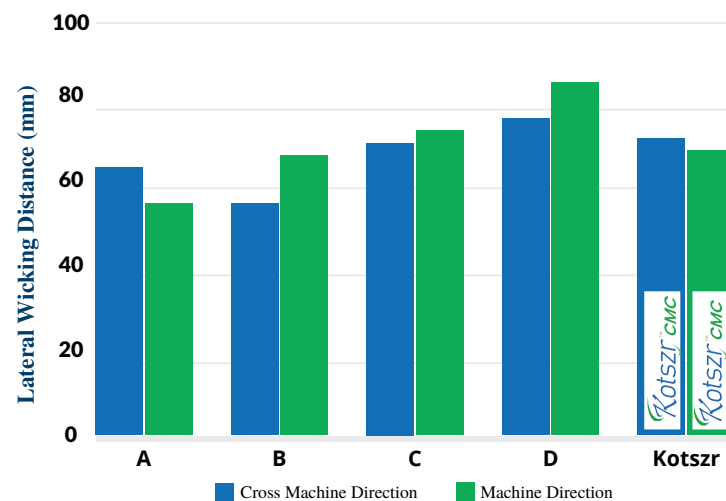

A - Aquacel is a registered trademark of ConvaTec

B - Aquacel Extra is a registered trademark of ConvaTec Inc.

C - Durafiber is a registered trademark of Smith & Nephew

D - Exufiber is a registered trademark of Mölnlycke

REFERENCES

1. L. Vinklářková, R. Mastejková, G. Folyňová, J. Muselík, S. Pavlovská, J. Bernatonič, D. Vetchý, Film wound dressing with local anesthetic based on in-soluble carboxymethylcellulose matrix, *J. Appl. Biomed.* 15 (2017) 313–320, <https://doi.org/10.1016/j.jab.2017.08.002>.
2. T. Maver, S. Hribnik, T. Mohan, D.M. Smrke, U. Maver, K. Stana-Kleinschek, Functional wound dressing materials with highly tunable drug release properties, *RSC Adv.* 5 (2015) 77873–77884, <https://doi.org/10.1039/C5RA11972C>.
3. M. Olsson, K. Järbrink, U. Divakar, R. Bajpai, Z. Upton, A. Schmidchen, J. Car, The humanistic and economic burden of chronic wounds: a systematic review, *Wound Repair Regen.* 27 (2019) 114–125, <https://doi.org/10.1111/wrr.12683>.
4. G.H. Kim, J.N. Im, T.H. Kim, G.D. Lee, J.H. Youk, S.J. Doh, Preparation and characterization of calcium carboxymethyl cellulose/chitosan blend nonwovens for hemostatic agents, *Text. Res. J.* 88 (2018) 1902–1911, <https://doi.org/10.1177/0040517517712101>.
5. P. Basu, U. Nanendrakumar, R. Anunachalam, S. Devi, I. Manjubala, Characterization and evaluation of carboxymethyl cellulose-based films for healing of full-thickness wounds in normal and diabetic rats, *ACS Omega* 3 (2018) 12622–12632, <https://doi.org/10.1021/acsomega.8b02015>.
6. M. Uzun, A review of wound management materials, *J. Text. Eng. Technol.* 4 (2018) 53–59, <https://doi.org/10.15406/jteft.2018.04.00121>.
7. G.D. Winter, Formation of the scab and the rate of epithelialization of superficial wounds in the skin of the young domestic pig, *Nature* 193 (1962) 293–294, <https://doi.org/10.1038/193293a0>.
8. I. Orlando, I. Roy, in: M.I.H. Mondal (Ed.), *Cellulose-based Hydrogels for Wound Healing*, Springer International Publishing, Cham 2018, pp. 1–18, https://doi.org/10.1007/978-3-319-76573-0_38-1.
9. A. Pettignano, A. Charlot, E. Fleury, Carboxyl-functionalized derivatives of carboxymethyl cellulose: towards advanced biomedical applications, *Polym. Rev.* 59 (2019) 510–560, <https://doi.org/10.1080/15583724.2019.1579226>.
10. T.W. Wong, N.A. Ramli, Carboxymethylcellulose film for bacterial wound infection control and healing, *Carbohydr. Polym.* 112 (2014) 367–375, <https://doi.org/10.1016/j.carbpol.2014.06.002>.
11. M.-C. Li, Q. Wu, K. Song, H.N. Cheng, S. Suzuki, T. Lei, Chitin nanofibers as reinforcing and antimicrobial agents in carboxymethyl cellulose films: influence of partial deacetylation, *ACS Sustain. Chem. Eng.* 4 (2016) 4385–4395, <https://doi.org/10.1021/acssuschemeng.6b00981>.
12. N.S.V. Capanema, A.A.P. Mansur, A.C. de Jesus, S.M. Carvalho, L.C. de Oliveira, H.S. Mansur, Superabsorbent crosslinked carboxymethyl cellulose-PEG hydrogels for potential wound dressing applications, *Int. J. Biol. Macromol.* 106 (2018) 1218–1234, <https://doi.org/10.1016/j.ijbiomac.2017.08.124>.

Why Choose Us?

- **PERFORMANCE**, significantly better in terms of absorbency.
- **STRENGTH**, helps in maintaining its structure from start to finish.
- **GENTLE**, facilitates intact removal and entire process trauma free.
- **PRICE**, with the aim to make advance wound care technologies accessible to all.

HOW TO APPLY DRESSING?

1



Always use aseptic techniques when dressing a wound

2



Clean and irrigate the wound

3



Open the *Kotszycmc* pouch & pull out the dressing. When folding or layering the dressing, make sure it fits the wound

4



Apply a secondary dressing and secure it

5



Irrigate with saline and gently peel off

6



Remove the dressing slowly and discard it

Special Characteristics

BIOCOMPATIBILITY


Unique de-watering process enhances the fibre strength, reduces fibre shed and provides increased biocompatibility.

ABSORBENCY

Optimal fibre spacing for increased absorbency. Conformable dressing that gels when in contact with wound exudate.

WET STRENGTH

Optimized needle bonding process designed to balance the conflicting requirements of strength and absorbency.

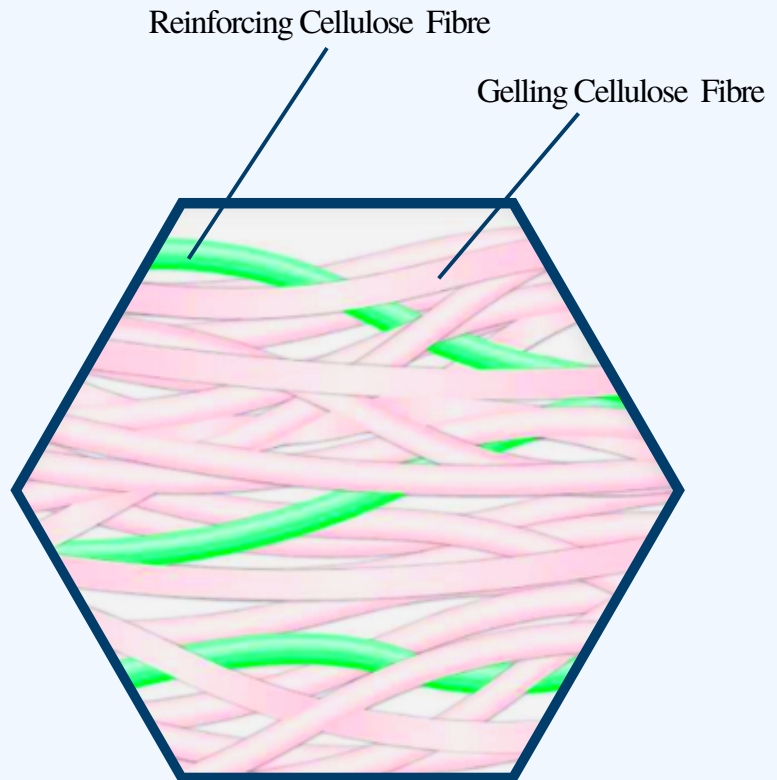

Ordering Information



Product Code	Size	Items
KZCMC0510	5×10 cm	10
KZCMC1010	10×10 cm	10
KZCMC1515	15×15 cm	10
KZCMC1020	10×20 cm	10
KZCMC1520	15×20 cm	10
KZCMC0230	2×30 cm	5
KZCMC3060	30×60 cm	5



Anziam Bio Private Limited

Manufacturing Facility: E-127(B), Bhiwadi Industrial Area, Near RIICO Chowk- Bhiwadi-301019 (Rajasthan)

R&D: Bionest Bio-incubator, Regional Centre for Biotechnology, Gurgaon- Faridabad Expressway, Faridabad-121001

91-7976625409

care@anziam.com

www.anziam.com

facebook.com/Anziam-Bio/

linkedin.com/company/anziam-bio/

https://bit.ly/Anziam-Bio